

Raucodrape®

Strong solutions, right down to the last detail - that is what the new product range is all about

Standards are defined by actual needs in daily practice

The strong points of Raucodrape are apparent in the smallest details - a multitude of innovative design features have been added for easier unfolding of the drape sheets, a wider variety of applications, and enhanced safety and convenience in surgical procedures.

Strong solutions in detail

■ Raucodrape® makes unfolding faster and easier

Raucodrape OP drapes now feature optimised folding styles. Depending on the specific purposes for which it is used, each type of drape has its own unique folding style which ensures safety in use and makes the process of unfolding faster and easier.

Now all our instrument table covers come with an extra-wide folding style for better handling and application.

In all self-adhesive fenestrated drapes, the aperture is visible even before the sheet is unfolded. This makes it considerably easier to apply the product correctly to the patient

■ Self-explanatory pictograms

Safe and fast unfolding of the drapes is also facilitated by the easily comprehensible pictograms. At a glance, they show the first step to unfolding the sheet, the direction in which it unfolds the correct positioning of the drape on the patient, as well as any splits, integrated leggings, and lots more.

■ Expanded range of single-sheet drape systems

The Raucodrape product range provides a comprehensive choice of modern single-sheet drape systems. The special features integrated in the single-sheet drapes save considerable time and effort in preparing for a surgical procedure.

As a general rule, all single-sheet drapes have a three-dimensional cut in the area of intervention, which prevents any undesirable tension and tautness in the drape. Thus, there is always enough space, even when several surgeons are working on a patient.

All integrated leggings in the single-sheet drapes are extra-long (125 cm), i.e. they can be used even with padded leg rests.

■ New integrated special features

L&R has placed increased focus on adding special features to


the Raucodrape OR drapes. The portfolio now includes even more of the following product categories:

- fluid collection pouches
- incision films
- tube retainers
- instrument pouches

■ Three-piece adhesive strips

Adhesive strips are an essential factor in the overall efficiency of any process, in particular with regard to product and user safety:

- All split sheet have an ergonomically designed, three-piece adhesive strip, which makes them easier to apply.
- Many adhesive sheets feature an extra-long adhesive strip to accommodate adipose patients.
- All adhesive strips are equipped with a 2 cm finger-lift for a better and safer grip.

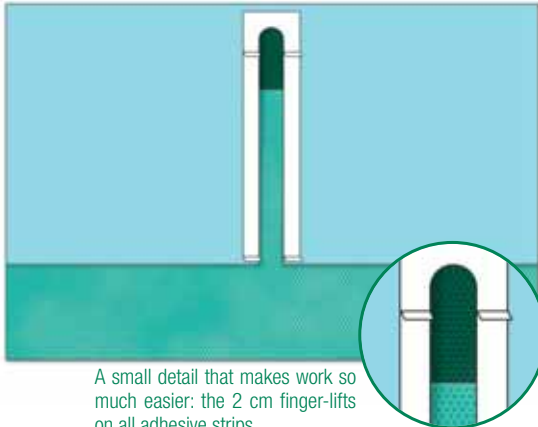

A small detail that makes work so much easier: the 2 cm finger-lifts on all adhesive strips